
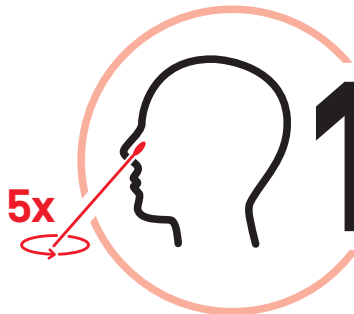


6 Step-by-step guide

on how to use the antigen self-test-kit

VIVADIAG PRO


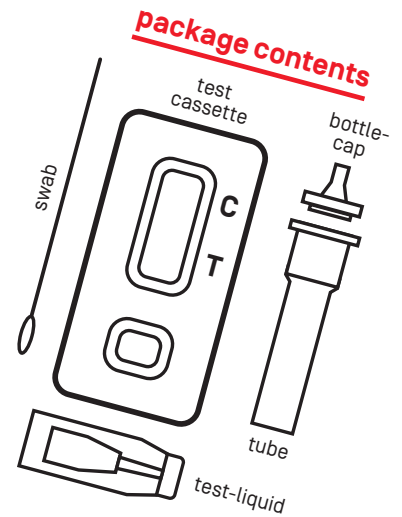


1.

Insert the swab **approx. 2,5 cm** into one nostril
Rotate the swab **5 times**, gently pressing it against the mucous membrane
Repeat the procedure in the other nostril

Hold the sealed capsule vertically and let the test liquid flow into the lower area
Break off the tip and pour all the liquid into the tube
Roll the swab **5 times**

2.

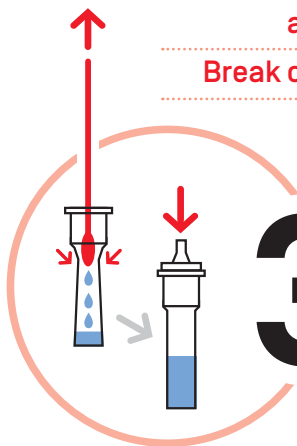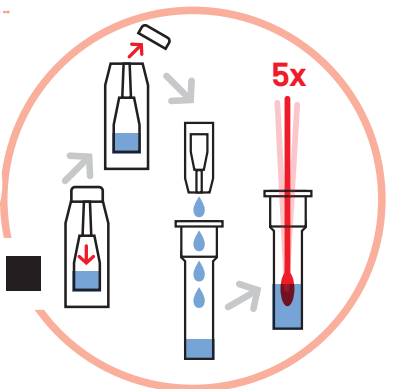


3.

Remove the swab from the tube and squeeze the tube at the same time
Press the cap with dropper top tightly on to the tube

Take the test cassette out of the packaging and place it flat on the table
Drop **3 drops** from the tube onto the test cassette

4.

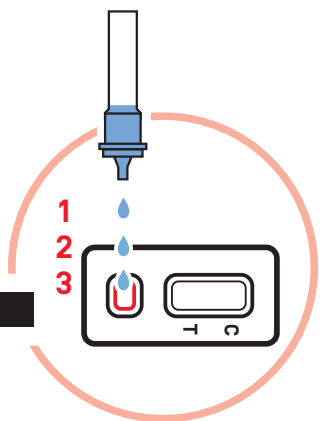

Important: The test cassette is sensitive to movement after the test liquid has been dripped onto it. In order not to falsify the test result, the cassette must therefore not be turned, shaken or, for example, carried from one room to another.




5.

Wait for **15 minutes**
Read the result

Interpret the result

6.

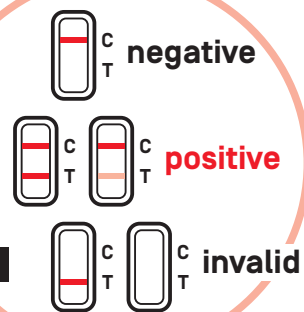

POSITIVE? → Call 1450!
NEGATIVE? Stay healthy!

Disclaimer: This guide is intended to provide a clear and easy-to-understand instruction on how to use the self-test kit to perform a rapid anterior-nasal COVID-19 antigen test. Please note that COVID-19 antigen rapid tests must always be carried out as specified in the instructions for use of the manufacturer or distributor. Please read the enclosed instructions for use carefully before using a COVID-19 antigen rapid test and follow the instructions. The Austrian Chamber of Pharmacists (Österreichische Apothekerkammer) assumes no liability for damage of any kind that may result from following the content of these instructions.

